


# Eclipse<sup>®</sup>

## Super Absorbent Dressing Range

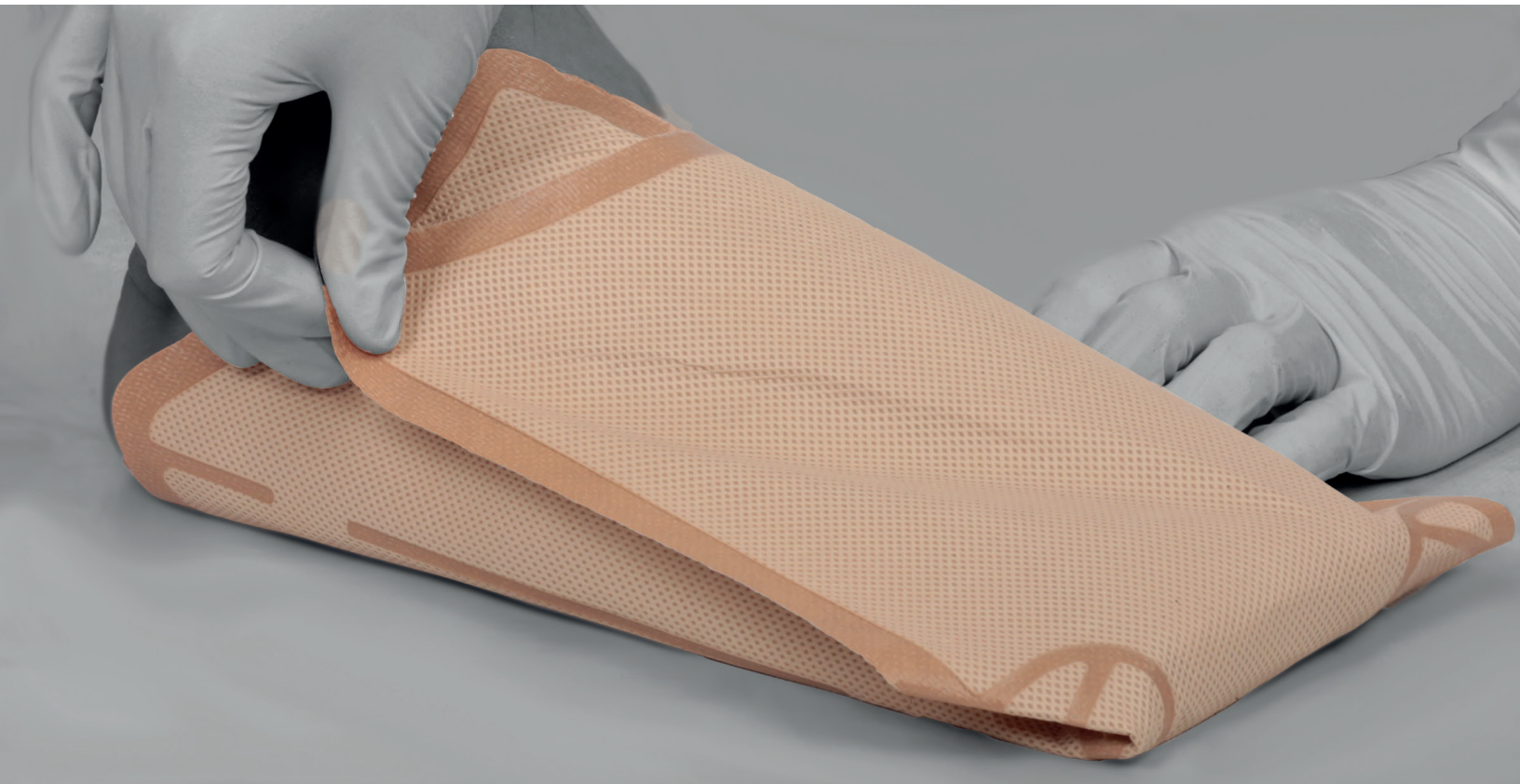

## Effective exudate management

Distributed by




[www.mediusa.com](http://www.mediusa.com)

**Adv**ancis<sup>®</sup>  
Medical

[www.advancismedical.com](http://www.advancismedical.com)

# How do the layers of the Eclipse® dressings work?

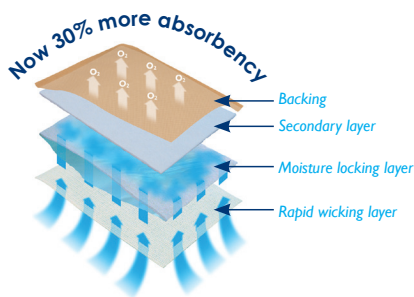

## The backing

The water-resistant barrier backing on Eclipse® dressings are bacteria and viral proof and critically prevent strike-through.



## Secondary layer

This layer sits underneath the backing and on top of the moisture locking layer.



## Moisture locking layer

The central layer of the Eclipse® dressing consists of a sheet of highly absorbent crystals and a mechanically bonded cellulose pad.



## Rapid wicking layer


This layer draws up the exudate from the wound bed delivering the exudate into the moisture lock layer.



## Soft silicone contact layer

The silicone contact layer gently adheres to skin but due to the hydrophobic nature of silicone it will not adhere to the wound bed.

Manufactured in our state-of-the-art facility in the UK, our range of Eclipse® super absorbent dressings deliver quality performance in the controlled management of exudate. Our recent developments in the Eclipse® range of dressings have seen an increase in the absorbency of our products by 30%, with no impact on cost.

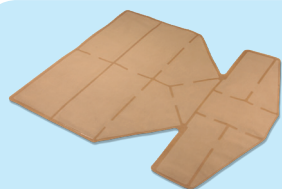

## Eclipse Boot®

Super absorbent leg wrap dressing

A unique clinical solution to assist in the management of excess fluid often associated with lower leg oedema.

Eclipse Boot® is a super absorbent dressing specifically designed to manage high levels of fluid in the lower limb.

The preshaped design overcomes the traditional time-wasting complexities of dressing 'leaky legs' with a combination of dressings that are often bulky and restrict patient mobility.

- ✓ Anatomically designed to fit the lower limb and foot
- ✓ 12 Super absorbent interconnected compartments to optimize absorbency
- ✓ Highly conformable
- ✓ Significantly improves patient mobility
- ✓ Easy to apply, single dressing design
- ✓ Reduce dressing application time
- ✓ Available in three sizes



## Eclipse Foot®

Super absorbent foot wrap dressing

A super absorbent dressing, designed to be applied securely around the foot area to manage wound exudate.

The dressing utilizes the high fluid handling capabilities of the Eclipse® range in a preshaped dressing designed to fit comfortably and effectively around the foot. Its is easy and simple to apply and offers enhanced patient comfort.

- ✓ Rapid absorption
- ✓ High capacity
- ✓ Fluid repellent backing
- ✓ Strength and durability
- ✓ Stay dry technology
- ✓ Thin and conforming
- ✓ Particularly useful in issues related to the treatment of the diabetic foot



Product	Size	Dressings per box	HCPCS code
Eclipse Boot® Small	22in x 19in	5	A6198
Eclipse Boot® Large	24in x 28in	5	A6198
Eclipse Foot® X-Large	28in x 31in	5	A6198
Eclipse Foot®	13in x 19in	5	A6198

Your patient outcomes can be improved when Eclipse® is used in conjunction with our range of 100% Medical Grade Manuka Honey products.



medi USA  
6481 Franz Warner Pkwy  
27377 Whitsett  
United States of America  
Tel: 800-633-6334 | E-mail: info@mediusa.com



www.advancismedical.com